

EXHIBIT G-2
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA SPECIALIST KELLY WOHLWEND AGENCY YAVAPAI CO SO

DATE 03/31/15 TIME 0704

INTOXILYZER SERIAL # 80-001095 LOCATION YCSO MAYER

() 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution, 100 AC.
OR

Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$

() 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"

() 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.

() 4. Air blank completed.

() 5. Calibration check completed. Test results 0. 097 AC.

() 6. Air blank completed.

() 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Kelly Wohlwend 9811

INTOXILYZER 8000
Location: YCSO-MAYER
Serial Number: 80-001095
Core Version: 8105.48
03/31/2015 02:04:51

Standard Lot#: 33913100A3
Last Changed By: L. MOODY #7486

QAS: K. WOHLWEND #9811
YAVAPAI CO SO

Operator: K. WOHLWEND #9811
YAVAPAI CO SO

Test	g/210L	Time
Air Blank	0.000	02:05:54
0.100 Cal Check	0.097	02:06:15
Air Blank	0.000	02:06:44

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-4
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST KELLY WOHLWEND AGENCY YAVAPAI CO, SO.
DATE 03/31/15 TIME 0208
INTOXILYZER SERIAL # 80-001095 LOCATION MAYER YCSO

1. Display Reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- 1. Clock time check.
- 2. Date check.

OPERATIONAL TESTS

- 1. Alcohol-free subject test result 0. 000 AC.
- 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- 4. Standard Calibration Check standard 0. 100 AC. Result 0. 096 AC.

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Kelly Wohlwend 9811

INTOXILYZER 8000

Location: YCSO-MAYER

Serial Number: 80-001095

Core Version: 8105.48

03/31/2015

02:08:56

Standard Lot#: 33913100A3

Last Changed By: L. MOODY

#7486

QAS: K. WOHLWEND

#9811

YAVAPAI CO SO

Operator: K. WOHLWEND

#9811

YAVAPAI CO SO

Subject: TEST,SQAP

DOB: 03/31/2000

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	02:09:56
Diagnostic Test	Pass	02:10:31
Air Blank	0.000	02:10:58
0.100 Cal Check	0.097	02:11:19
Air Blank	0.000	02:11:49
Subject Test	0.000	02:12:27
Air Blank	ABT*	02:12:34
Air Blank	IPS**	02:12:41

*Sequence Aborted

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: YCSO-MAYER

Serial Number: 80-001095

Core Version: 8105.48

03/31/2015

02:13:57

Standard Lot#: 33913100A3

Last Changed By: L. MOODY #7486

QAS: K. WOHLWEND #9811

YAVAPAI CO SO

Operator: K. WOHLWEND #9811

YAVAPAI CO SO

Subject: TEST,SQUAP

DOB: 03/31/2000

Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	02:14:59
Diagnostic Test	Pass	02:15:34
Air Blank	0.000	02:16:01
0.100 Cal Check	0.097	02:16:22
Air Blank	0.000	02:16:51
Subject Test	0.023*	02:19:59
Air Blank	0.000	02:20:29
Air Blank	0.000	02:20:56
Subject Test	XXX**	02:21:16
Air Blank	RFI***	02:21:27

*Deficient Sample

**Invalid Sample

Begin new deprivation period

***RFI Detect

Not a Successfully Completed Test Sequence

EXHIBIT G-2
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

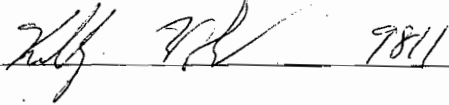
STANDARD CALIBRATION CHECK PROCEDURE

QA SPECIALIST KELLY WOHLVEND AGENCY YAVAPAI CO. SO.

DATE 03/31/15 TIME 0221

INTOXILYZER SERIAL # 80-001095 LOCATION YCSO MAYER

- () 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution .100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- () 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- () 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- () 4. Air blank completed.
- () 5. Calibration check completed. Test results 0.096 AC.
- () 6. Air blank completed.
- () 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

INTOXILYZER 8000

Location: YCSO-MAYER

Serial Number: 80-001095

Core Version: 8105.48

03/31/2015

02:21:52

Standard Lot#: 33913100A3

Last Changed By: L. MOODY #7486

QAS: K. WOHLWEND #9811

YAVAPAI CO SO

Operator: K. WOHLWEND #9811

YAVAPAI CO SO

Test g/210L Time

Air Blank 0.000 02:22:23
0.100 Cal Check 0.096 02:22:44
Air Blank 0.000 02:23:13

Preventative Maintenance
Performed -- Timer Reset