

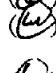






**EXHIBIT G-2  
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)**

**ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000**

**STANDARD CALIBRATION CHECK PROCEDURE**

QA SPECIALIST Logan Moody #7486 AGENCY Yavapai County Sheriff's Office  
 DATE 12/19/2018 TIME 1259  
 INTOXILYZER SERIAL # 80-000289 LOCATION YCSO - Black Canyon

-  1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 100 AC.  
 OR  
 Pour a standard alcohol concentration solution ~~AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made.~~ Turn on the simulator and allow ~~temperature to reach 34°C ± 0.2°C~~
-  2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
-  3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
-  4. Air blank completed.
-  5. Calibration check completed. Test results 0. 099 AC.
-  6. Air blank completed.
-  7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 