



INTOXILYZER 8000

Location: PRESCOTT P.D

Serial Number: 80-000337

Core Version: 8105.48

11/15/2016

03:56:34

-----  
Standard Lot#: 22315100A1

Last Changed By: S. HILTON #372  
-----

QAS: S. HILTON #372

PRESCOTT PD  
-----

Operator: S. HILTON #372

PRESCOTT PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 03:57:08 |
| 0.100 Cal Check | 0.099  | 03:57:29 |
| Air Blank       | 0.000  | 03:57:59 |

Preventative Maintenance  
Performed -- Timer Reset

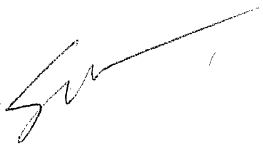


EXHIBIT G-4  
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST S. Hertz AGENCY Prescott P.D.  
DATE 11/15/16 TIME 0357  
INTOXILYZER SERIAL # 80-000737 LOCATION 222 S. MARION ST.

6/5 1. Display Reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- 6/5 1. Clock time check.  
6/5 2. Date check.

OPERATIONAL TESTS

- 6/5 1. Alcohol-free subject test result 0. 000 AC.  
6/5 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
6/5 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
6/5 4. Standard Calibration Check standard 0. 100 AC. Result 0. 090 AC.

Instrument is operating properly and accurately. YES  NO

COMMENTS

\* FUNCTIONING PROPERLY \*

SIGNATURE S. Hertz

INTOXILYZER 8000  
Location: PRESCOTT P.D  
Serial Number: 80-000337  
Core Version: 8105.48  
11/15/2016

03:59:23

---

Standard Lot#: 22315100A1  
Last Changed By: S. HILTON #372

---

QAS: S. HILTON #372  
PRESCOTT PD

---

Operator: S. HILTON #372  
PRESCOTT PD

---

Subject: SQAP, TEST  
DOB: 12/19/1973  
Sex: M Weight: 175

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 04:00:37 |
| Diagnostic Test | Pass   | 04:01:11 |
| Air Blank       | 0.000  | 04:01:38 |
| 0.100 Cal Check | 0.098  | 04:01:59 |
| Air Blank       | 0.000  | 04:02:29 |
| Subject Test    | 0.000  | 04:02:45 |
| Air Blank       | 0.000  | 04:03:12 |
| Wait            | ABT*   |          |
| Air Blank       | 0.000  | 04:03:47 |

\*Sequence Aborted

*SMA*

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: PRESCOTT P.D

Serial Number: 80-000337

Core Version: 8105.48

11/15/2016

04:04:23

---

Standard Lot#: 22315100A1

Last Changed By: S. HILTON #372

---

QAS: S. HILTON #372

PRESCOTT PD

---

Operator: S. HILTON #372

PRESCOTT PD

---

Subject: SQAP, TEST

DOB: 12/19/1973

Sex: M Weight: 175


---

15 Minute Deprivation Period? Yes

---

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | IPS*   | 04:05:08 |
| Air Blank | 0.000  | 04:05:37 |

\*Improper Sample

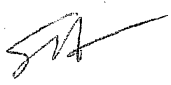

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: PRESCOTT P.D  
Serial Number: 80-000337  
Core Version: 8105.48  
11/15/2016

04:06:01

---

Standard Lot#: 22315100A1  
Last Changed By: S. HILTON #372

---

QAS: S. HILTON #372  
PRESCOTT PD

---

Operator: S. HILTON #372  
PRESCOTT PD

---

Subject: SQAP, TEST  
DOB: 12/19/1973  
Sex: M Weight: 175

---

15 Minute Deprivation Period? Yes

---

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | RFI*   | 04:06:43 |
| Air Blank | 0.000  | 04:07:11 |

\*RFI Detect

*SQA*

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: PRESCOTT P.D  
Serial Number: 80-000337  
Core Version: 8105.48  
11/15/2016

04:07:32

---

Standard Lot#: 22315100A1  
Last Changed By: S. HILTON #372

---

QAS: S. HILTON #372  
PRESCOTT PD

---

Operator: S. HILTON #372  
PRESCOTT PD

---

Subject: SQAP, TEST  
DOB: 12/19/1973  
Sex: M Weight: 175

---

15 Minute Deprivation Period? Yes


---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 04:08:45 |
| Diagnostic Test | Pass   | 04:09:19 |
| Air Blank       | 0.000  | 04:09:47 |
| 0.100 Cal Check | 0.099  | 04:10:08 |
| Air Blank       | 0.000  | 04:10:38 |
| Subject Test    | 0.204* | 04:13:45 |
| Air Blank       | 0.000  | 04:14:30 |
| Air Blank       | 0.000  | 04:14:57 |
| Subject Test    | XXX**  | 04:15:20 |
| Air Blank       | 0.000  | 04:16:23 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

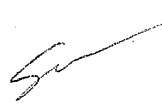

**Not a Successfully Completed Test Sequence**

EXHIBIT G-2  
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA SPECIALIST S. HEITW AGENCY PROSECUTOR P.D.

DATE 11/15/16 TIME 0916

INTOXILYZER SERIAL # 80-000337 LOCATION 222 S. MARINA ST.

(S) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 1.00 AC.  
OR

Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$

- (S) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"  
(S) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.  
(S) 4. Air blank completed.  
(S) 5. Calibration check completed. Test results 0. 0.90 AC.  
(S) 6. Air blank completed.  
(S) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE S. HEITW

INTOXILYZER 8000

Location: PRESCOTT P.D

Serial Number: 80-000337

Core Version: 8105.48

11/15/2016

04:16:56

-----  
Standard Lot#: 22315100A1

Last Changed By: S. HILTON #372  
-----

QAS: S. HILTON #372

PRESCOTT PD  
-----

Operator: S. HILTON #372

PRESCOTT PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 04:17:28 |
| 0.100 Cal Check | 0.090  | 04:17:49 |
| Air Blank       | 0.000  | 04:18:19 |

Preventative Maintenance  
Performed -- Timer Reset

